

L I G H T I N G



MIDMARK

Where many lights are an afterthought of old technology, the Midmark light pushes lighting quality to the edge. It features an integrated PC board that monitors and adjusts the voltage to ensure color corrected brightness, LED intensity readouts and touchpad controls. All of this helps to deliver the true, reduced shadow white light you need to properly match shades, identify details and diagnose tissue. The feathered edge of the light pattern minimizes eye strain while the durable polycarbonate shield protects the reflector and is easily removed for cleaning and bulb replacement. User-friendly features simplify using the light while minimizing fatigue, time and motion.




The touchpad control offers three intensity adjustments as well as a composite safe setting to prevent premature curing. The automatic motion controlled on/off function can be set on the touchpad with the simple touch of a button.




The Track Light is centrally located and offers an extended range of motion for greater flexibility. An optional HD monitor and flexarm allow for seamless integration of PC, TV and DVD.



The Ceiling Mounted Light houses the transformer under the mount cover, providing maximum operator clearance with a modern, uncluttered design.

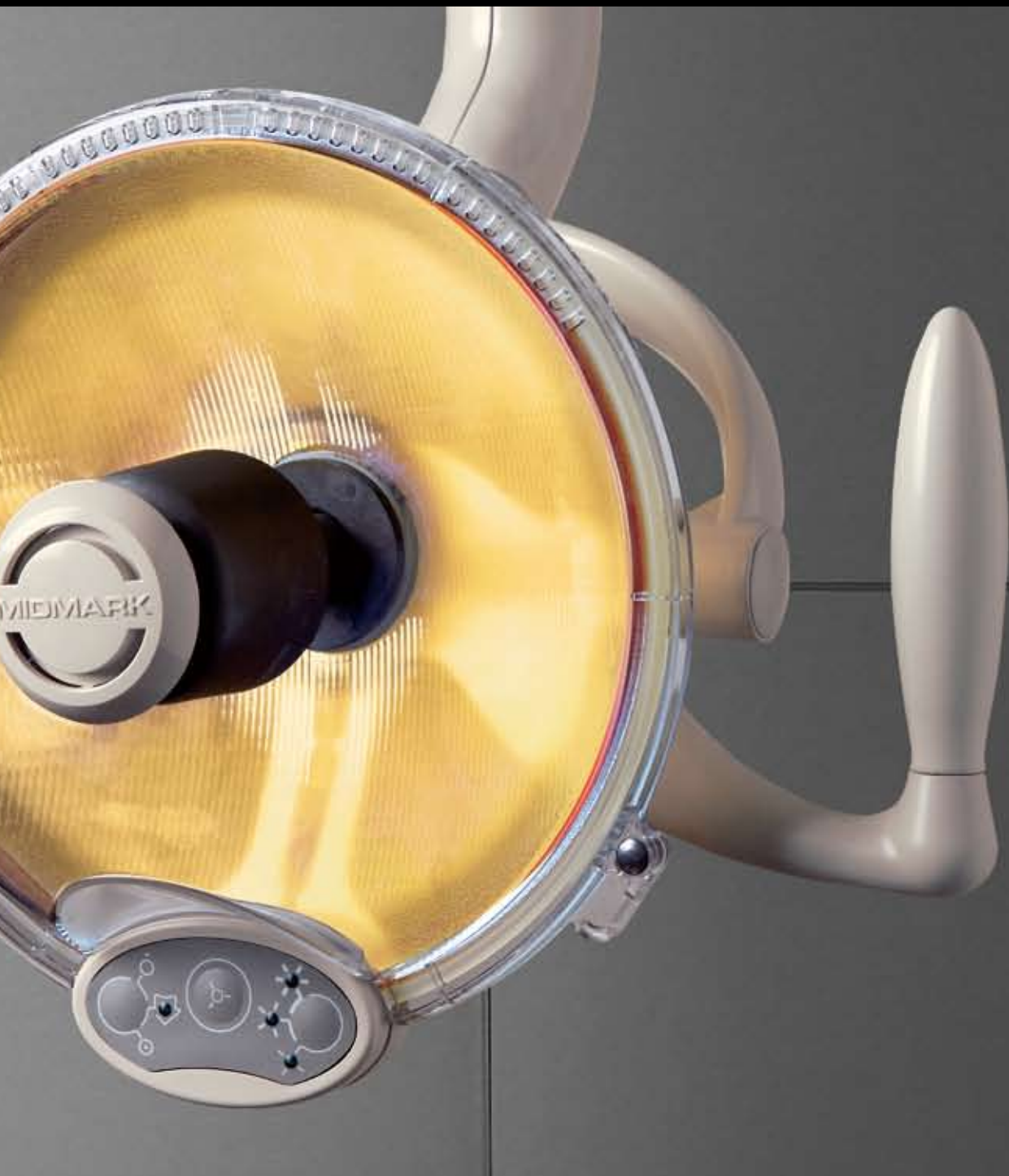

The Wall and Cabinet Mounted Lights offer excellent stability and allow the light to be conveniently positioned against the wall or cabinet when not in use.

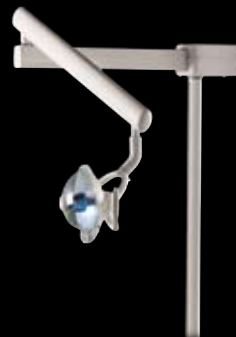

The Universal Mounted Light houses the transformer on the light arm to provide maximum adaptability for chair mounting. This light will fit round posts with bearings inside diameters of 1 3/4" or 1 3/8".

Lighting Specifications

Operate Mode Color Temperature:	4500, 4750, 5000 Kelvin
Operate Mode Intensity:	1400, 1950, 2500 Foot Candles 15000, 21000, 27000 Lux
Composite Mode Color Temperature:	4000 Kelvin
Composite Mode Intensity:	550 Foot Candles 5900 Lux
Color Rendering Index:	Greater than 90%
Pattern Size:	3.3" x 7.5"

Midmark is an ISO 13485 and ISO 9001 Certified Company.
For more information or a demonstration, contact your Midmark
dealer or Call: 1-800-MIDMARK Fax: 1-877-725-6495
Outside the U.S.A. Call: 1-937-526-3662 Fax: 1-937-526-8214
or visit our website at midmark.com

© 2009 Midmark Corporation
60 Vista Drive P.O. Box 286 Versailles, Ohio 45380-0286
Products subject to improvement changes without notice
Litho in U.S.A. 007-0720-00 Rev. A1 (8/09)



Because we care.